

BioLogic.

Your one-stop-shop for accessories and add-ons.



BioLogic Add-on instruments

Rotating Electrodes: BluRev	Systems Tips Cells Enclosure	06-09
Temperature control units	Intermediate Temperature System (ITS-e) Controlled Environment Sample Holder (CESH-e) Conductivity cell (HTCC) High Temperature Furnace (HTF-1100) High Temperature Sample Holders (HTSH-1100)	10-13

Accessories

Battery accessories	Coin Cell Holders Cylindrical Cell Holders Pouch Cell Holders Prismatic & Pouch Cell Holders Current Collectors Sense Adapter Module (SAM-50) Redox Flow Battery cells (RFB)	16-21
Analytical cells	Small volume cells Large volume cells Multi purpose cells	22-25
Corrosion cells	Standard corrosion cells Avesta cell Flat cells Galvanic cells Plate material evaluating cell Coating cells	26-29
Electrodes	Working electrodes Counter electrodes Small reference electrodes Hydrogen reference electrodes	30-35
Connection accessories	High Temperature Cables Glove box cables Multi-electrode investigation cables Connectors Connection kits External device connection Test boxes Faraday cage	36-38
Spectroelectrochemistry	Spectrometer system Spectroelectrochemical Cell	39-41


Your partner from A to Z. Turnkey solutions provider



Master the **full measurement chain** using high quality instruments and accessories. The quality, reliability and innovation that lies behind BioLogic's product portfolio helps shape the future of research and industry around the world. Take advantage of the **unique specs** of your instruments.

Master the full measurement chain

Our range of add-on instruments and accessories are designed to increase the scope of your electrochemical experiments, without compromise: corrosion, material study, electrochemical reactions and much more can be easily and fully addressed.

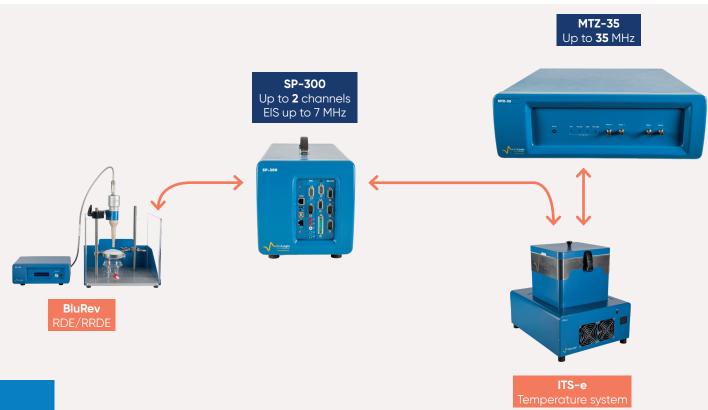


Battery: from research to industry

Throughout the full battery value chain every single component must be thoroughly tested. This includes anode, cathode, binder, separator, electrolyte all the way to the commercial cell, and there are unique challenges at every step:

- Materials & components research
- Research cell performance & characterization
- Manufacturing process optimization
- · Commercial cell validation & ageing
- Screening & benchmarking for modules
- Integration & second life evaluation



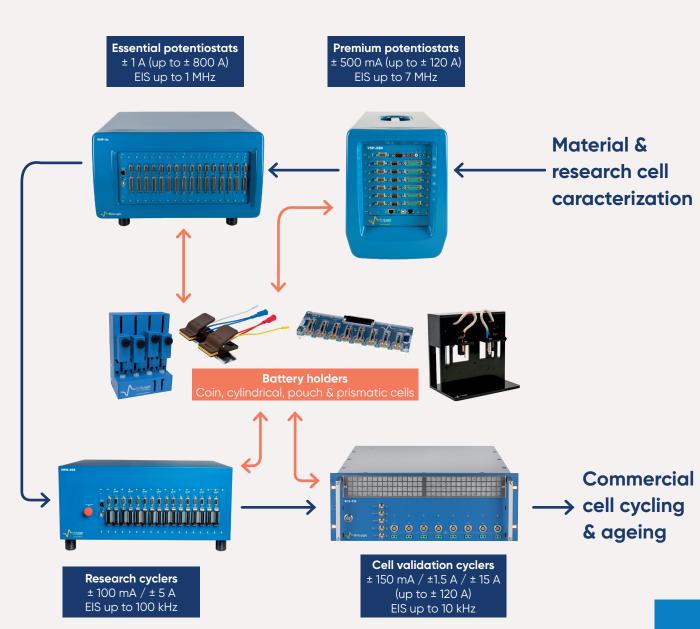

BluRev Rotating disk/ring disk electrodes

Steady state efficiency



Exploit the power of EC-Lab® with BluRev rotating electrodes

BluRev is a range of robust, versatile rotating disk (RDE) and rotating ring disk (RRDE) electrodes ideal for use with BioLogic potentiostats. A range of quick-fit exchangeable electrode tips (multiple materials/diameters) can be easily added to the body. Driven by EC-Lab® software, an industry standard for potentiostat control software, BluRev instruments can be used independently, or as part of an integrated setup.

A fully integrated solution

The BluRev's custom designed enclosure has been specifically designed to facilitate the operation and set-up of the BluRev RRDE. This makes it easy to hold the RRDE body as well as the BioLogic potentiostat cables. The result is a fully integrated, purpose-built rotating electrode system.





BluRev Systems

±1 RPM constant accuracy over the whole rotational range

These modular research instruments demonstrate excellent levels of accuracy, particularly at high rotation rates. The RC-10k control unit offers an accuracy of ±1 rpm over the whole rotational range for precise and fully reproducible experimental conditions.



Access the true rotation rate even in very demanding media

BluRev systems are the only devices on the market to display both target and true rates of rotation thanks to the embed optical encoder. Accessing this data can greatly assist users during experiments by enabling them to monitor rotation rates in real-time. Viscous medias will not affect the rotation speed as the system will adapt it accordingly.

Robust, even in highly corrosive conditions

Protect the rotating electrode by connecting it to inert gases. The gas will flow through the shaft inner parts to protect them from corrosion.

Highly compatible

The speed can be set manually or remotely by using the analog output of a BioLogic instrument.




094-RDE without tip

Specifications

	Information
Rotational range/rpm	100 - 10,000 (9000 for RRDE)
Accuracy/rpm	1 typical over the whole rotational range
Runout/mm	<0.1
Setting resolution/rpm	10 (Manual control) or 25 (Remote control)
Materials of RRDE/RDE	PEEK, Al alloy
Inert gas inlet for shaft corrosion protection/mm	Ø 2
Operating temperature/°C	10 - 40
Power	24 Vdc, 1 A max
Max consumption/W	24
Dimensions/mm	RC-10k: 95x227x178 (HxWxD), RDE/RRDE: 233.6 (length with tip)
Weight/kg	RC-10k: 1.00 RDE/RRDE: 0.36 (without tip)

	094-RC/RDE	094-RC/RRDE
Content		
RC-10k Rotation controller	09	4-RC
Electrode rotator (motor, shaft, electrode body, Ag/C brushes)	094-RDE	094-RRDE
DB9 to BNC connector for external control of RC-10k	092-22/1	
1 m BNC/BNC cable	Included	
Replacement Ag/C brush	094-RDE-BRUSH	
1 transport case	Included	

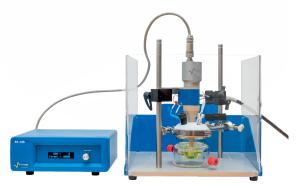

BluRev Tips

All standard tip bodies (M6 thread) are made of PEEK, but for experiments requiring a high chemical resistance, PTFE versions are also available.


All tips are polished to obtain a final roughness Ra of 50 nm (*Ra = 10 nm for the Boron Doped Diamond tip).

RDE Tips

	Catalog n°
Glassy Carbon (3 mm)	094-GC/3
Glassy Carbon (5 mm)	094-GC/5
Glassy Carbon, PTFE body (3 mm)	094-PTFE-GC/3
Glassy Carbon, PTFE body (5 mm)	094-PTFE-GC/5
Platinum (2 mm), 99.9 %	094-Pt/2
Platinum, PTFE body (3 mm), 99.9 %	094-PTFE-Pt/3
Gold (2 mm), 99.9 %	094-Au/2
Gold, PTFE body (3 mm), 99.9 %	094-PTFE-Au/3
Boron Doped Diamond (3 mm), 500-1000 ppm 500 µm diamond layer	094-BDD/3
Titanium (3 mm), 99.9 %	094-Ti/3
Silver (3mm), 99.9 %	094-Ag/3
Aluminum (3 mm), 99.9 %	094-AI/3
Copper (3 mm), 99.9 %	094-Cu/3
Nickel (3 mm), 99.9 %	094-Ni/3
Stainless steel (3 mm), 99.9 %	094-316L/3

RRDE Tips

		Catalog n°
Glassy C (Nmax =	arbon Ring & Disk 0.4626)	094-GC-GC
Pt Ring (9 (Nmax =	9,9 %) - Glassy Carbon Disk 0.4626)	094-Pt-GC

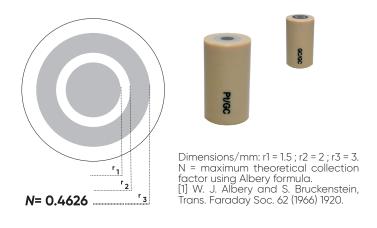



Replaceable Disk Tips

This solution allows users to put custom designed materials as working electrodes (disks) to perform RDE experiments.

A PEEK empty tip is available to perform experiments on 3 mm OD disks. The material is placed inside the empty tip by using a mounting tool kit and spacers: leak free solution.

	Catalog n°
Mounting tool kit	094-DRE-KIT
3 mm PEEK empty tip	094-DRE/3
Spacer kit (1x PEEK, 3x PTFE)	094-SPACER-KIT
Polishing tip	094-DRE-POLISH

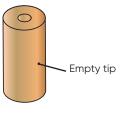

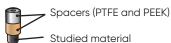

BluRev Cells

The BluRev RDE/RRDE range is compatible with a wide range of cells, in particular the

EL-ELECTRO-80, EL-ELECTRO-80DJ or EL-ELECTRO-150DJ kits.

	Catalog n°
PTFE 5 holes cap compatible with BluRev RDE (needed if you already have an EL-ELECTRO cell)	094-A-CAP

Note: if you already possess one of these cells, you will need a special adaptor cap: **094-A-CAP**.

BluRev Enclosure



	Catalog n°	
Protective housing and stand kit for the BluRev	094-ENCL	
Contents		
1x stainless steel plate with M6 threaded holes to fix support poles 2x support poles to hold BluRev RDE and the cell 1x clamping flange for the BluRev 3x half-clamps for the potentiostat cables (all BioLogic cables are supported) 1x plexiglas protective housing		
Specifications		
Dimensions with protective housing (HxWxD)/mm	287x318x308	
Weight (with protective housing)/kg	5.3	

BluRev Background information and theory

RDE


Commonly used techniques...

Levich and Koutecký-Levich are powerful analysis tools used to obtain kinetic electrochemical parameters such as the diffusion coefficient of a redox species in a given medium and the reaction constant (Application Note #56)

...advanced EIS based techniques

Fitting impedance measurements made on a redox reaction occurring at a rotating electrode at only one rotation speed also enables the direct measurement of the diffusion coefficient.

By Using ZFit and W_{inf} (BioLogic unique element), it is possible to directly obtain the diffusion coefficient of the species of interest.

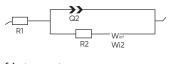

For more detailed information please see the EC-Lab application note #66.

-Im(Z) vs. Re(Z) — RDE_Z_Fell_Felll_SmM_01_PEIS_zfit.mpp • RDE_Z_Fell_Felll_SmM_01_PEIS.mpr, cycle 3

Nyquist plot

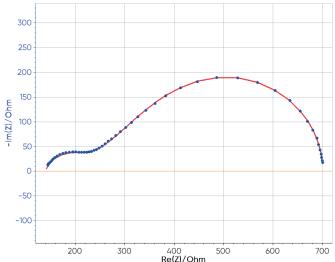


Figure 1: PEIS measurement using an equimolar solution of 5 mM $\rm K_3Fe(CN)_6$ and $\rm K_4Fe(CN)_6$ in 0.1 M KCl and a BluRev RDE at 2000 rom.

RRDE

BioLogic bipotentiostats

At the disk electrode, the electroactive species are generated depending on the applied potential and detected at the ring. In a typical experiment, a CV is performed on the disk electrode and a constant voltage is applied on the ring electrode. EC-Lab® software embed a CV-CA technique that is dedicated for this purpose ("bipotentiostat" techniques folder).

For RRDE measurements, a bipotentiostat is needed to control both working electrodes. Because of the presence of two working electrodes, a specific connection mode is needed to avoid ground loops: the "CE to Ground" connection mode.

Unique "CE to Ground" connection mode for accurate measurements

The "CE to Ground" connection mode is recommended as it highly reduces risks of current leakage that leads to measurement errors. All BioLogic multichannel potentiostats offer this feature.

Temperature Control Unit.

Master sample temperature



ITS-e

A small-footprint benchtop temperature chamber

The intermediate temperature system allows accurate temperature control between -35°C and +150°C. This is a good alternative to the larger conventional climatic chambers available on the market. Thanks to its small chamber, the temperature homogeneity is much better than in conventional chambers and leads to more reproducible measurements.

Application oriented

The ITS-e is dedicated to control the material's temperature for its electrical characterization. It can also be used for thermal curing and aging and for the determination of some kinetic or thermodynamic parameters like activation energy, phase transformation or glass transition temperatures.

Fully integrated and compatible with BioLogic products

The ITS-e is compatible with the MTZ-35 impedance analyzer and with BioLogic potentiostats. It allows users to easily set up their electrical and electrochemical experiments thanks to EC-Lab® or MT-Lab® software.





Specifications

	Information
Temperature range /°C	-35 to +150
Temperature accuracy /°C	± 0.3
Input voltage	115 V / 230 V, 50/60 Hz
Computer Interface	USB or Ethernet
Power consumption /W	250
Dimension /mm	400 x 313 x 385 (HxWxD)
Weight /kg	8
Catalog n°	
ITS-e	097-140e

CESH-e

A leak-tight sample holder for controlled atmosphere study

The Controlled Environment Sample Holder is a leak-tight (up to 2 bar relative) sample holder for the electrical characterization of hard, powdered, pasty and soft materials. Oxygen and moisture sensitive samples can be prepared in a glove box and then placed inside the CESH-e.

The CESH-e can be used as standalone, or inside an ITS-e or a climate chamber.

Constant contact pressure and thickness measurement

A thickness measurement kit including a micrometer head with a calibrated ratchet is provided as an option to ensure the accurate measurement of the thickness of the sample (ASTM D374). This also enables reproducible pressure (100 N or 6 N respectively) to study flexible and compressible materials.



	Information
Operating temperature /°C	-40 to 150
Max sample diameter /mm	30
Max sample thickness /mm	4
Core material	Anodized aluminum
Electrode material	Gold plated copper
Dimension /mm D 79 x H 94 mm	79x94 (Diameter x Height) 121.90 clearance
Residual capacitance with 20 mm diameter /pF	8
Weight /kg	0.9
To be used	Standalone or with ITS-e or climate chamber
Catalog n°	
CESH-e	097-150e
Thickness measurement kit	097-150e/10
Cable kit for CESH-e and MTZ-35	097-150e/01
Cable kit for CESH-e and Potentiostat	097-150e/02

Electrodes

Thanks to its modular electrodes setup, the CESH-e allows Through-plane and In-plane electrical measurements using interchangeable electrodes.



(Left) In-plane / (Right) Trough-plane

Trough plane
097-150/TP20 (20 mm diameter) 097-150/TP12 (1/2" diameter) 097-150/TP06 (1/4" diameter)

HTCC

High Temperature Conductivity cell (up to 180°C)

The HTCC is a sealed 2-pole conductivity cell composed of two parallel platinum electrodes. Its design allows measurements from -50°C up to 180°C. With a nominal constant of 1.0 cm⁻¹, conductivities from $2\,\mu\text{S}\cdot\text{cm}^{-1}$ to $0.2\,\text{S}\cdot\text{cm}^{-1}$ can be measured.



Platinized and non-platinized electrodes

The platinized electrodes are coated with a black platinum layer in order to increase the effective surface area of the electrode. Only 0.5 mL of liquid/gel is required, reducing the cost of experiments and their environmental impact. The cells are compatible with the ITS-e.

	Information
Cell type	Platinum parallel plates on glass holder
Connections	2-wire (I = 70 cm each)
Nominal constant /cm-1	1 ± 10 %
Minimal sample volume /mL	0.5
Conductivity range /µS.cm	2 to 200 000
Temperature range /°C	-50 to +180
Maximum temperature ramp for platinized cell /°C/min	1 (to avoid platinum detachment)
Compatible with	Potentiostat, MTZ-35, ITS-e
Catalog n°	
Platinized (x1)	098-010/10
Non-platinized (x1)	098-010/11

HTF-1100

High Temperature Furnace (up to 1100°C)

The HTF-1100 is a horizontal laboratory tube furnace dedicated to the electrical characterization of materials and for heat treatment in the range between 150°C and 1100°C.

Manual and remote control

The HTF-1100 can be controlled by using the Watlow controller or MT-Lab® software. Its controller facilitates the set-up and consequent monitoring of the temperature during tests. The controlled temperature can be based on a k-type thermocouple placed in the bottom of the furnace or one placed inside the sample holder (HTSH-1100): accurate control and measurement of the temperature close to the sample.



	Information
Temperature range /°C	150°C up to 1100°C
Temperature control accuracy /°C	< ± 1
Temperature ramp /°C/min	Adjustable from 0.1 to 20
Insulation material	Alumina Fiber
Temperature sensors	K-Type Thermocouple
Safety Features	Emergency stop button Buzzer sound alarm Temperature safety limit
Catalog n°	
HTF-1100	097-111

HTSH-1100

The furnace accommodates multiple tubular High Temperature Sample Holders (HTSH-1100). It can operate under controlled environment conditions with inert or active gas (Ar, N_2 , O_2 , etc) and with variable pressures up to 2 bar relative thanks to its gas inlet/outlet and quartz tube and safety valve.

	Information
Temperature range /°C	150°C up to 1100
Quartz tube dimension /mm	45.5x270.5 (Diameter x Height)
Electrodes diameters /mm	3 / 12 / 25
Electrodes material	High purity Platinum disc (1 mm thickness)
Maximum sample diameter /mm	27
Maximum sample thickness /mm	5
Temperature probe	K-Type thermocouple (inconel shielded)
Water cooling	Tube included for water cooling when using T > 800°C
Dimension (including quartz tube) /mm	80x310.5 (Diameter x Height)
Weight /kg	1.2
Catalog n°	
HSTF-1100, 3 mm diameter electrode	097-130/S
HSTF-1100, 12 mm diameter electrode	097-132/S
HSTF-1100, 25 mm diameter electrode	097-133/S


TCU Background information and theory

Battery component researches


Researching components for batteries often requires precise data over a large temperature range to determine how temperature affects ionic conductivity. With precise electric conductivity measurements collected over a range of temperatures, EC-Lab® is able to calculate the material's activation energy, and thus determine the temperature stability of a material.

Advanced temperature control

A solid electrolyte sample is loaded inside a CESH-e and put inside an ITS-e to collect multiple impedance measurements at various temperatures. EC-Lab® software is used to control both VMP-300 potentiostat and ITS-e by using the TCU server.

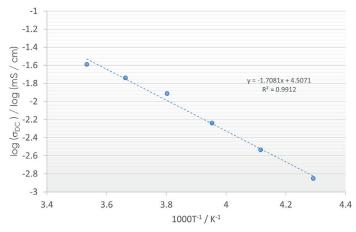

Arrhenius plot for the bulk conductivity of solid electrode: Temperature intervals of 10°C from -40°C to 10°C

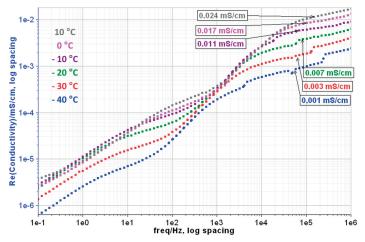

AC conductivity obtained from EIS data: Temperature intervals of 10° C from -40° C to 10° C

Determine the activation energy (Arrhenius plot)

By setting the sample's thickness and diameter, EC-Lab® can calculate the sample cell constant (k) and directly use it to plot the electrical conductivity as a function of frequency: DC conductivities can then be extracted from the pseudo-plateau found in the high frequency domain (shown above). By using this information, the activation energy for the materials conductivity can be determined from an Arrhenius plot, shown on the left (0.14 eV).

A low activation energy indicates that the material's ionic conductivity will remain stable when exposed to temperature extreams and the material can be used for applications in a wider range of temperature environments.

Premium.

When only the best will do



VMP-300 16 channels, ± 150 A down to 100 fA



SP-200 1 Channel, ± 500 mA down to 100 fA



SP-300 2 Channels, ± 10 A down to 100 fA



VSP-300 6 Channels, ± 40 A down to 100 fA

BioLogic Premium is a range of state-of-the-art potentiostats designed for researchers who need the fastest, most precise, potentiostats available.

Premium potentiostats boast some of the most powerful specifications available: 100 fA to 150 A, 7 MHz EIS measurements and a sampling rate that can reach $1 \text{ data point every } \mu \text{s}$.

Premium by name. Premium by nature.

Essential.

Tools for **Electrochemists**

VMP-3e, 16 channels, ± 1 A down to 20 nA





SP-50e 1 channel, ± 1 A down to 20 nA



SP-150e 2 channels, ± 1 A down to 20 nA



VSP-3e 8 channels, ± 1 A down to 20 nA

A range of powerful, modular, high-precision, research-grade potentiostats built to handle almost any academic or industrial application imaginable.

From 1 to 16 channels. 1 A (native), up to 800 A with boosters and EIS QI™.

Essential measurement tools. Whatever your area of specialisation.




Battery Accessories. Four point cell holders Ensure reliable, repeatable measurements

Four points are better than two

By measuring only the impedance of the cell, a 4-point connection battery holder enables reliable and repeatable measurements. Its design negates measurements of connector and holder-related impedances.

All of our battery holders, with the exception of the CCH model, are built around the 4-point connection design.

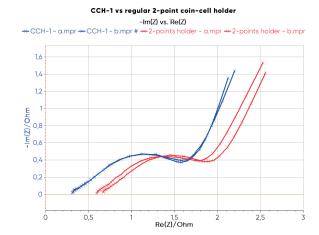

Coin Cell Holders







	CCH-1	CCH-8	ССН		
Cell max diameter/mm		24			
Cell height/mm		1.6 - 3.2	3		
Number of channels	1	8	4		
Measurement type 4/2 point		2 point			
			MPG2		
To be used with	A	ny instruments VMP3, VMP-3e			
Climatic chamber compatibility	Yes (-30 to 80°C)	Yes (-30 to 80°C)	No		
Catalog n°	096-126	096-128/H Cables must be purchased separately	092-22/14		

CCH-8 Cables



Compatible with	Length	Catalog n°
BCS-905	2.5 m	096-128/C25
BCS-910	5 m	096-128/C50
Essential/Premium range	75 cm	092-22/24a
Potentiostats MPG2 (banana receptacles)	2.5 m	092-22/24b
Essential range Potentiostats MPG2 (front panel connexion)	2.5 m	092-22/24c



Cylindrical Cell Holders



BH-1i

4-point measurement

CBH-8

4-point measurement

2 and 4 mm receptacles are available for the current (power) cables. For voltage (sense) cables, only 2 mm receptacles are available.

4 mm receptacles are available for the current (power) and voltage (sense) cables.



CBH-4

	BH-1i	CBH-4	СВН-8	
Cell max diameter/mm	26	60		
Cell min height/mm	0	3	0	
Cell max height/mm	76	10	00	
Number of channels	4	4	8	
Max current/A	12	3	52	
Measurement type	4 point	4 p	oint	
Receptacles diameter/mm	2 and 4	2	4	
To be used with		All instrument	S	
Max operating T/°C	60	8	0	
Size : HxWxD/mm	205x150x95	335x260x150	335x520x150	
Weight/kg	0.6	1.9 3.8		
Catalog n°	092-22/15	092-C32/4 092-C32/8		

More reliable measurements, higher repeatability levels.

Pouch Cell Holders



PBH-125

	PBH-125	PBH-150	PBH-4	PBH-8	
Min leads separation distance/mm		0	12		
Max leads separation distance/mm		110*	4	44	
Number of channels		1	4	8	
Max current/A	25	50	-	32	
Measurement type		4 point			
Receptacles diameter/mm	4 (power) 2 (voltage)	6 (power)** 4 (voltage)	4		
To be used with		All instrume	nts		
Max operating T/°C	80	100	80		
Size : HxWxD/mm	40x5	50×210***	135x325x180 135x650>		
Weight/kg	0.2*** 1.9		3.8		
Catalog n°	092-P25/1	092-P50/1	092-P32/4	092-P32/8	

*Measured using the guide rail and the middle of the clamp.

PPRH-1100

^{**}Eyelet ring (connection kit **094-110/CNT** can be used for an easier connection to 6 mm diameter cables).

***Measured with the two clamps mounted on the guide rail.



PBH-4

4-point measurement

measurement

PBH-4 & PBH-8 holders

4 mm receptacles are available for the current (power) and voltage (sense) cables.

Prismatic and Pouch Cell Holders



PPBH-1100

4-point measurement

	PPDN-132	PPBH-1100		
Cell min height/mm	0			
Cell max height/mm	139)		
Min leads separation distance/mm	66			
Max leads separation distance/mm	155	5		
Number of channels	1			
Max current/A	32	100		
Measurement type	4 point			
Receptacles diameter/mm	4 (power and sense) 4 (power and sense) (power)			
To be used with	All instru	ments		
Max operating T/°C	80)		
Size : HxWxD/mm	265x320x300 320x320x360			
Weight/kg	3 5.1			
Catalog n°	092-PC32/1	092-PC100/1		

PPRH-132

Pouch cell and prismatic holder

4 mm receptacles can be used for currents up to 32 A. For higher currents, the 6 mm flush mounting plugs should be used. These are compatible with FlexP 0160, HCV-3048, CC4-60 A and CC8 cables.



Pouch cell and prismatic holder

4 mm receptacles are available for the current (power) and voltage (sense) cables.

Current Collectors

BioLogic's current collectors offer the possibility to connect several channels in parallel and increase the maximum current that can be passed through the cell, in order to simplify and reduce the footprint of your setup.







	CC4-60A	CC8	CC4-200A		
Connection details					
Input					
Power cables/receptacles diameter/mm	4	/ 4	6 (IP2x)		
Voltage sense receptacles diameter/mm	2	2	4		
Number of input channels	4	8	4		
Max current/channel/A	1	5	50		
Output					
Power receptacles diameter/mm	6 (1	P2x)	8 (Amphenol, IP2x)		
Voltage sense receptacles diameter/mm	2 (16	P2x)	4 (IP2x)		
Max output current/A	60	120	200		
Cables details					
Output power cables	1 pair of 2 m power cable	es with 6 mm receptacles	1 pair of 2.5 m power cables with 8 mm receptacles and M8 threads		
Output voltage cables	1 pair of 2 m sense cables	with 2 mm banana plugs	1 pair of 2.5 m sense cables with 4 mm banana plugs		
Instrument compatibility	VSP-	-915* -300 -300	FlexP0160 FlexP0060 HCV-3048		
Included connection kit	094-11	0/CNT	093-200/CNT		
Measurement type		4 p	point		
Max operating T/°C		80			
Size (with feet) : HxWxD /mm	70x170x88	70x300x88	120×248×169		
Weight/kg	3.8				
Catalog n°	096-022	096-015	093-100/CC4		

^{*}The CC8 comes with BCS tablets and cables. It is also compatible with the VSP-300 and VMP-300 and can be provided without cables and tablets using the following part number: **096-015/1.**

Sense Adapter Module (SAM-50)



This can be added to a multi-channel system to perform stack measurements up to 60 V for 5 channel boards and a 10-element measurement. 3 SAM-50s can be linked to track up to 30 elements.

Sense Adapter Module	Catalog n°
SAM-50	092-26

Redox Flow Battery Cells






BioLogic provides two types of redox flow batteries with two different active surface areas: 5 cm² and 20 cm². This package allows you to work with both aqueous or organic electrolytes depending on your research needs or studies. Ready to use out of the box, we have different testing packages available that consist of 8 felts and 4 membranes. You can choose between Fumasep F-1850, mutliple Nafion types or VANADion membranes.


Sealings and bipolar plates can also be purchased separately as spare parts.

Complete turnkey setups with climatic chambers and flow control are also available. Please ask your local reseller for more information or see the table below.

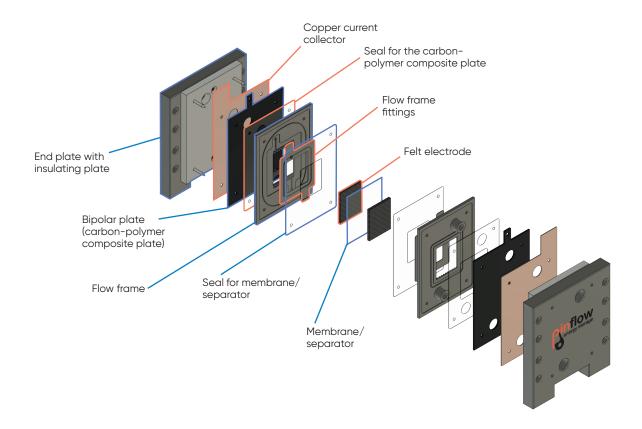


Pinflow RFB cells exhibit low internal resistance as can be seen on the following curve. Typical values of internal resistance are lower than $2 \Omega \cdot \text{cm}^2$.

Pinflow cells were used and characterized in the following papers: https://doi.org/10.1016/j.memsci.2018.02.011 https://doi.org/10.1016/j.jpowsour.2018.01.079

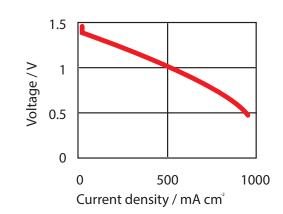

The Lab-cells from Pinflow are specially designed to control the pressure applied on the carbon felts that are used as electrodes. Using rigid components and non-flat bipolar plates, one can not only perform reproducible experiments, but also, it is possible to use various electrode thicknesses.

The schematic below shows, that by providing non-flat carbon polymer composite plates, one can control the compression and the thickness of the carbon felt electrode being used. The thickness of the electrode compartment can be easily set up using carbon-polymer composite plates with a defined stump or pit.

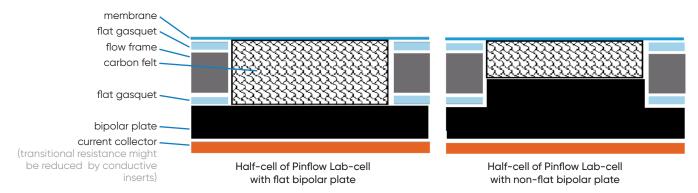

			5 cm ²	20 cm ²	
Lab Cells		Standard	P-LABCELL/5	P-LABCELL/20	
		For organic media	P-LABCELL/5F	P-LABCELL/20F	
Consumables	Testing packages	Fumasep F-1850	P-TEST-PACK/5-11	P-TEST-PACK/20-11	
	(8x felts, 4x membranes)	Nafion N115	P-TEST-PACK/5-13	P-TEST-PACK/20-13	
		Nafion N117	P-TEST-PACK/5-14	P-TEST-PACK/20-14	
		Nafion N212	P-TEST-PACK/5-15	P-TEST-PACK/20-15	
		Vanadion	P-TEST-PACK/5-16	P-TEST-PACK/20-16	
	Membranes (x1)	Fumasep F-1850	P-MEMB/5-11	P-MEMB/20-11	
		Nafion N115	P-MEMB/5-13	P-MEMB/20-13	
		Nafion N117	P-MEMB/5-14	P-MEMB/20-14	
		Nafion N212	P-MEMB/5-15	P-MEMB/20-15	
		Vanadion	P-MEMB/5-16	P-MEMB/20-16	
Spare parts	Pack of sealings (x10)	Standard	P-PSEAL/5	P-PSEAL/20	
		For organic media	P-PSEAL/5F	P-PSEAL/20F	
	Sealings	Standard	P-SSEAL/5	P-SSEAL/20	
		For organic media	P-SSEAL/5F	P-SSEAL/20F	
	Bipolar plates (1x set of 2)	-	P-SBIPOL/5	P-SBIPOL/20	
	Tubing and vessels	-	P-TUBINGS		
	Fittings	-	P-FPM		
Accessories	Felt cutter	-	P-CUTTER/5	P-CUTTER/20	
	Wrench	-	P-WRENCH		
	Stand	-	P-STAND		
	Peristaltic pump	-	P-LABCELLPUMP/2		
	Climate chamber	-	P-CCHAMBER		
	Data acquisition module	_	P-DATAMOD		

Application Note: Need More Information?

In this application note, a Vanadium Redox Flow Battery (VRFB) was characterized using typical DC and AC techniques: galvanostatic charge and discharge cycling and Electrochemical Impedance Spectroscopy (EIS).



Analytical Cells.

Small Volume Cells

Each voltammetry cell is designed for a specific application (specific working electrode, volume of solution, oxygen-free condition, etc.). As an example, for voltammetry investigations using standard working electrodes with an outer diameter (OD) of 6 mm, the fixed configuration of the SVC-3 kit is recommended. For applications requiring other working electrode shapes, the SVC-2 is more suitable.

6 mm, the fixed configuration of the Si recommended. For applications requiorking electrode shapes, the SVC-2 is m




The SVC-2 in microvolume mode was designed for those cases where only a small amount of electroactive compound is available.

Here is a list of the available cells:

- SVC-2, modular
- SVC-3, for a volume of 5 to 20 mL, only for working electrode with OD of 6 mm
- VC-4, for a volume of 1 to 3 mL, only for working electrode with OD of 6 mm
- Bulk electrolysis cell, for a volume of 100 mL




							1	
		Catalog n°		Catalog n°		Catalog n°		Catalog n°
Products	SVC-2	A-012668	SVC-3	A-012669	VC-4	A-011224	Bulk electrolysis cell	A-001197
Content								
Sample vial/mL	20 (7 pieces)	A-001056	20 (7 pieces)	A-001056	5 (7 pieces)	A-011504	100 (1 piece)	A-012632
Counter electrode (CE)/ mm	57	A-002222	50	A-002233	57	A-002222	230	A-002234
PTFE cap		A-012670		A-012671		A-011226		A-012551
Purge tube (ETFE), 100 mm	l	Included		Included		Included		Included
Additional items	Adapter 10 to 6 mm	Included			Cell holder	A-011227	Porous carbon elec- trode	A-010530
							Lid for CE	A-001198
							Chamber for CE	A-001196
							O-ring	A-001236
							Port plug	A-009131
							Stirrer bar	A-000178
Options								
Sample holder/mm	9.0 (2 pieces)	A-012177	6.0 (2 pieces)	A-012176				
Cell holder	for 20 mL	A-001209	for 20 mL	A-001209				
Purge tube (ETFE)/m	1	A-010537	1	A-010537	1	A-010537	1	A-010537
Working electrodes				(See page 26			
Reference electrodes	See page 28							
PTFE cap	10.2 mm 1.6 mm 6.2 mm 1.6 mm 1.6 mm							

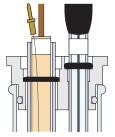


Small amount of solution (200 μ L in the sample holder)


SVC-2 modularity: microvolume mode

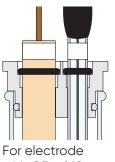

with OD of 10 mm

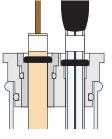

For electrode with OD of 6 mm SVC-2 modularity: oxygen free mode

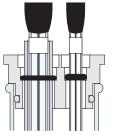

with OD of 4 mm

30 20 10 -10 /nA -20 -30 -40 -50 -0.5

Please note that a full, purpose-built, analytical kit is also available SK-2 (A-012763) but the reference electrode must be purchased separately. This kit includes:

- SVC-3 kit (A-012669)
- PK-3 polishing kit (A-011975) see page 26
- one glassy carbon electrode, OD 6.0 mm, ID 3.0 mm (A-002012) see page 26
- one platinum electrode, OD 6.0 mm, ID 1.6 mm (A-002013) see page 26

Cell Geometry


 E_{we}/V

The geometry of the cell should be optimized to reduce the ohmic drop. Working and reference electrodes must be positioned close to one another. The counter electrode should not limit the transfer of electrons, so its contact surface should be larger than the contact surface of the working electrode.



Large Volume Cells

This cell is perfect for standard analytical electrochemistry experiments with electrolyte volumes of several tens of mL. It is also compatible with the BluRev product range, *i.e.* electrode rotators.

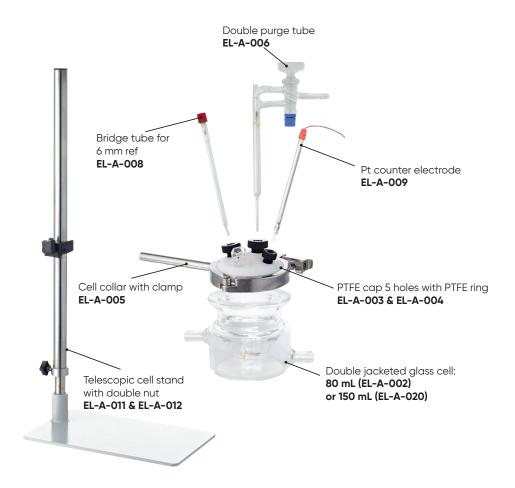

	EL-ELECTRO-80	EL-ELECTRO-80DJ	EL-ELECTRO-150DJ	
Glass cell	EL-A-001 (80 mL)	EL-A-002 (80 mL double jacketed)	EL-A-020 (150 mL double jacketed)	
PTFE cap 5 holes		EL-A-003		
PTFE ring, silicon encapsulated, OD 102 mm		EL-A-004		
Cell collar with clamp		EL-A-005		
Bridge tube for reference electrode, OD 6 mm		EL-A-008		
Platinum counter electrode	EL-A-009			
Purge tube	EL-A-016	-		
Double purge tube	- EL-A-006			
Double nut 25 mm and 12 mm diameter	- EL-A-011			
Telescopic cell stand	- EL-A-012			
Options				
Electrode bridge extension for electroanalytical cell	EL-A-022			
Bridge tube for reference electrode of OD 8 mm		EL-A-017		
PT100 probe	EL-C-014			
220 V - Magnetic stirrer & header, without PT100	EL-C-015A			
110 V - Magnetic stirrer & header, without PT100	EL-C-015B			
Aluminum base holder for magnetic stirrer	EL-C-018			
Set of 10 porous 4 mm glass frits (CoralPor) with PTFE heat shrink (200 mm)		092-VYC4		

Multi Purpose Cells

FlexCell®

These cells, manufactured by Gaskatel, are ideal for corrosion experiments in aggressive media, as well as studies on Gas Diffusion Electrodes (GDE) and membranes. Their unique design and a specific choice of materials allow repeatability and avoid common pitfalls in other cells: heterogeneous electric field, variable ohmic drop, crevice corrosion, degradation of the cell. Used in combination with the robust HydroFlex® or MiniHydroFlex® hydrogen reference electrode, this is the perfect cell for corrosion, membrane and GDE studies. Two versions are available: PTFE and PP.





FlexCell® PTFE

FlexCell® PP

	Catalog n°
FlexCell® PP – Electrochemical Test Cell made of PP	G-FLEXCELL/PP
FlexCell® Analyte Compartment in PP for membrane studies	G-COMP/PP
FlexCell® PTFE – Electrochemical Test Cell made of PTFE	G-FLEXCELL/PTFE
FlexCell® Analyte Compartment in PTFE for membrane studies	G-COMP/PTFE
Options to be ordered separately (see p.31)	
HydroFlex® Hydrogen Reference Electrode	G-HYDROFLEX
HydroFlex® Hydrogen Reference Electrode Starter Kit	G-HYDROFLEX-KIT
MiniHydroFlex® Hydrogen Reference Electrode	G-MINIHYDROFLEX






PTFE analyte compartment

PP analyte compartment

	Operating conditions
Max. exposure time/h	24
Temperature/°C	-20 to 120
pH range	-2 to 16
Max. current/A	3

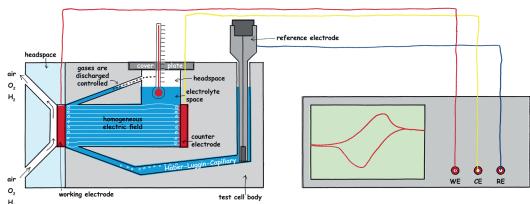

	Max. sample size /cm	Max. sample thickness /mm*	Active area /cm²	Electrolyte volume /ml	Dimensions (w/o screw adapters) (HxWxD)/mm
FlexCell®	3 x 5	10	3	30	75 x 100 x 132
Analyte Compartment	3 x 5	10	3	12	75 x 100 x 30

Cross-section of the cell and instrument connection.

^{*}Longer wing screws can be supplied

	Materials
Main	PP or PTFE
Gasket	Silicon
O-rings	EPDM
Cover plate	PSU
Gas compartments	PMMA or PSU
Counter electrode	Pt-Ir


PP	Polypropylene
PTFE	Polytetrafluoroethylene
PSU	Polysulfone
PMMA	Polymethylmethacrylate
EPDM	Ethylene propylene diene monomer
Pt	Platinum
lr	Iridium

Corrosion Cells.

Standard Corrosion Cells

The standard corrosion cell is available with as single or double jacket glass cell for temperature control. It is provided with two graphite rods to be used as counter electrodes, a bridge tube, to ensure the minimum distance between the reference and the working electrode, and purge tubes, to maintain a controlled gaseous environment.

Two kits are available: a standard one and an advanced kit with telescopic cell stand, sample holder and double purge tubes.

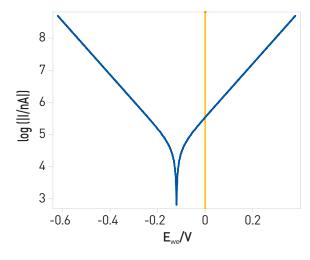




		Basic corrosion cell kit EL-CORR-1	Complete corrosion cell kit EL-CORR-1DJ	
Glass cell 1 L		EL-C-001	EL-C-002 (double jacketed)	
PTFE cap		EL-C	- 003	
PTFE ring, silicon encapsulated, OD 102 mm		EL-C	-004	
Cell collar with clamp		EL-C	:-005	
Graphite counter electrode rod (2 pieces) ρ = 1.070 $\mu\Omega$ cm		EL-C	-009	
Bridge tube for 6 mm diameter reference electrod	е	EL-C-008		
Purge tube		EL-C-016	-	
Double purge tube		-	EL-C-006	
Double nut 25 mm and 12 mm diameter		-	EL-C-011	
Telescopic cell stand		-	EL-C-012	
Sample holder 1 cm 2 (max sample thickness 3.4 mr 14.6 mm)	m and max diameter	- 092-C-016		
Options				
Bridge tube for 8 mm diameter reference electrod	е	EL-C-017		
PT100 probe (indicate connector type)		EL-C-014		
Magnetic stirrer & heater, without PT100 probe 220 V 110 V		EL-C-015A		
		EL-C-015B		
Aluminum base holder for magnetic stirrer and 1 L	cell vial	EL-C-018		
Set of 10 porous frits (4 mm CoralPor™) with PTFE h	neat shrink (200 mm)	092-	VYC4	

Avesta Cell


The Avesta Cell is an electrochemical cell developed for pitting corrosion testing (ASTM G150). It is designed to avoid microcrevice corrosion formed between the working electrode and the gasket at the bottom aperture of the cell.

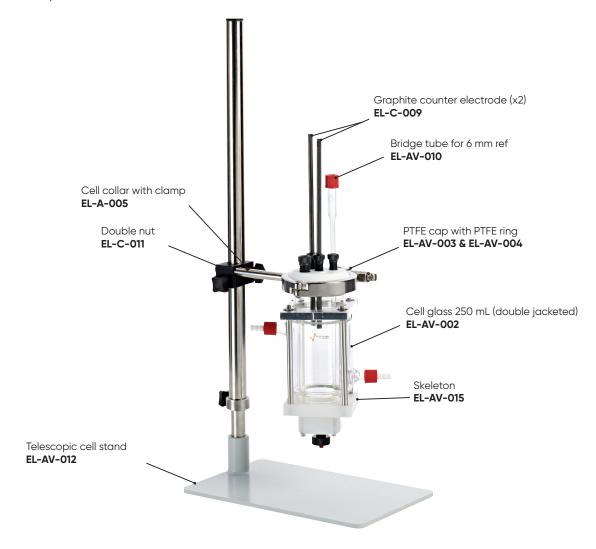

	Catalog n°
Avesta cell kit	EL-AVESTA
Content	
Double jacketed cell glass 250 mL	EL-AV-002
PTFE cap 5 holes	EL-AV-003
O-ring PTFE silicone encapsulated	EL-AV-004
Cell collar with clamp	EL-A-005
Double purge tube	EL-AV-006
Filter paper ring (100 pieces)	EL-AV-007
Graphite counter electrode rods (2 pieces) ρ = 1.070 $\mu\Omega$ cm	EL-C-009
Bridge tube for RE 6 mm	EL-AV-010
Double nut	EL-C-011
Telescopic cell stand	EL-AV-012
Skeleton	EL-AV-015
Options	
Peristaltic pump for low flow	EL-AV-008
Bridge tube for reference electrode with OD of 8 mm	EL-AV-013
Single purge tube	EL-AV-014
Temperature probe PT100	EL-C-014
Set of 10 porous frits (4 mm CoralPor™) with PTFE heat shrink (200 mm)	092-VYC4

A filter paper ring placed between the sample and the gasket is flooded by distilled water in order to eliminate crevice corrosion.


The water flow is controlled by a peristaltic pump (EL-AV-008).

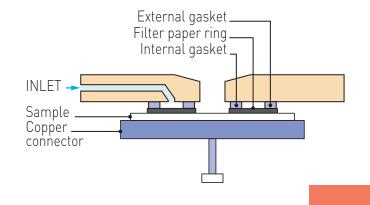

Flat Cells 1 to 10 cm² sample area

The Flat cell with a volume of 250 mL is perfect for experiments on flat specimens of 1 or 10 cm² surface area.

This cell has a double jacket for temperature control and three holes for reference electrodes, purgetubes and temperature probes with an inner diameter of 17.6 mm and the two others of 8.3 mm.

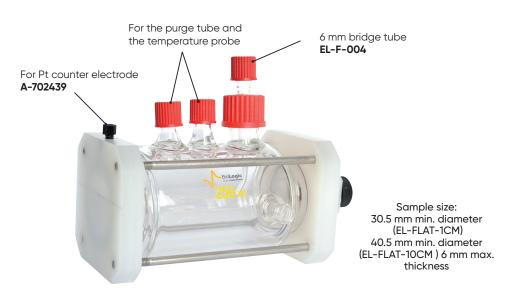

		s s kits e jacket)	PVDF kits (single jacket)		
	1 cm ² EL-FLAT-1CM	10 cm ² EL-FLAT-10CM	1 cm ² EL-FLAT-1CM-PVDF	10 cm ² EL-FLAT-10CM-PVDF	
Platinum mesh counter electrode (54 mm wire/80 mesh), 25 x 35 mm	A-702439				
Bridge tube for 6 mm diameter reference electrode	EL-F	-004		-	
Glass cell (250 mL) for flat cell	EL-F	EL-F-002			
Mechanical parts	EL-FLAT-3H EL-FLAT-4H		EL-FLAT-3H	EL-FLAT-4H	
Option					
Bridge tube for 8 mm diameter reference electrode	EL-F-	-004B		-	

Investigations In Agressive Media

If the experiment is performed in more agressive media such as hydrofluoric acid, it is possible to obtain the body of the flat cell in PVDF* instead of glass (Polyvinylidenefluoride).

	-	ofluoric d 48%		furic d 98%		phoric d 85%		chloric d 35%		itric d 70%		hloric cid		dium kide 50%	hyd	roxide acentr.
Temperature	20° C	50° C	20° C	50° C	20° C	50° C	20° C	50° C	20° C	50° C	20° C	50° C	20° C	50° C	20° C	50° C
PTFE																
PVDF*																
Borosilicate glass																
														exce		
Options									Ca	talog n°				not		ended
PVDF body (single-	jackete	d)						EL-F-P	VDF						compat	

Options	Catalog n°
PVDF body (single-jacketed)	EL-F-PVDF
EPDM O-rings for 1 cm ² with PEEK ferrule for CE*	EL-SEAL-1B
EPDM O-rings for 10 cm ² with PEEK ferrule for CE*	EL-SEAL-10B
PTFE O-rings for 1 cm ² with PEEK ferrule for CE*	EL-SEAL-T1B
PTFE O-rings for 10 cm ² with PEEK ferrule for CE*	EL-SEAL-T10B
EPDM O-rings for 1 cm ^{2*}	EL-SEAL-1
EPDM O-rings for 10 cm ^{2*}	EL-SEAL-10
PTFE O-rings for 1 cm ^{2*}	EL-SEAL-T1
PTFE O-rings for 10 cm ^{2*}	EL-SEAL-T10
Set of 10 porous frits (4 mm CoralPor™) with PTFE heat shrink (200 mm)	092-VYC4

 $^{^{\}star}$ The O-ring kits include 4 O-rings for the glassware side and 10 O-rings for the sample side

Galvanic Cells 1 to 10 cm² sample area

Thanks to the modular design of the flat cell, it is possible to place two different materials at each end of the cell. The surface area can be 1 or 10 cm². This cell is provided with a double jacket by default.

	1 cm ² EL-GAL-1CM	10 cm ² EL-GAL-10CM
Content		
Flat cell kit 1 cm ²	EL-FLAT-1CM	EL-FLAT-10CM
Galvanic kit 1 cm ²	092-FLAT/1	092-FLAT/10



Sample size:
30.5 mm min. diameter (EL-GAL-1CM)
40.5 mm min. diameter (EL-GAL-10CM)
6 mm max. of thickness

Plate Material Evaluating Cell up to 1 cm² sample area

This cell was developed to evaluate plate material such as metals, semi-conducting plates,

The sample plate is sandwiched between the two cell blocks. The required volume of solution is about 1 mL.



	Catalog n°
Plate material evaluating cell	A-011951
Content	
PTFE cell [body & base] (1 piece)	Included
O-ring (1 piece)	Included
Screw 20 mm (2 piece)	Included
Purging tube, 100 mm	Included
Platinum counter electrode (1 piece)	A-002222
Options	
O-ring (10 pieces)	A-012022

Coating Cell

The "Coating cell" is an affordable cell especially dedicated for testing flat coated material samples. Masks with different areas are available for corrosion testing.



	Catalog n°
Coating cell kit	EL-COAT
Content	
Glass for coating cell	EL-P-002
Nylon base with three feet	EL-P-003
Rubber cup with two holes	EL-P-004
Metallic clamp	EL-P-005
O-ring for coating cell	EL-P-006
Graphite rod counter electrode (L: 145 mm, OD: 6 mm, ρ = 1.070 $\mu\Omega$ cm)	EL-P-009
Options	
Bridge tube for 6 mm reference electrode	EL-P-008
Mask for 1 cm ² (20 pieces)	EL-P-011
Mask for 3 cm² (20 pieces)	EL-P-012
Mask for 10 cm ² (20 pieces)	EL-P-013

Electrodes.


Working Electrodes

To address every application, a wide range of working electrodes (WE) is available with diameters ranging from 7 μ m up to 10 mm and made of different materials.

BioLogic exclusive



M-BDD-3: Boron-doped diamond 3 mm diameter disk

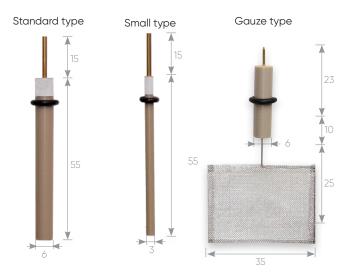

Dimensions in mm


			Isolation	OD/mm	Disk Ø/mm	Catalog n
Carbon	Standard glassy carbon electrode		PEEK	10	5	A-002417
			PEEK	6	3	A-002012
		PEEK	6	1.6	A-012297	
	Small glassy carbon electrode		PEEK	3	1.6	A-012298
			PEEK	3	1	A-002412
	Standard pyrolytic graphite electrode	Basal plane	PEEK	6	3	A-002252
		Edge plane	PEEK	6	3	A-002253
	Standard plastic formed carbon electron	ode	PEEK	6	3	A-002408
			PEEK	6	1	A-002409
	Small plastic formed carbon electrode		PEEK	3	1	A-011854
Platinum (99.95% purity)	Standard platinum electrode		PEEK	10	5	A-002420
			PEEK	6	3	A-002422
		PEEK	6	1.6	A-002013	
Small platinum electrode			PEEK	3	1.6	A-002313
Gold	Standard gold electrode		PEEK	10	5	A-002418
			PEEK	6	3	A-002421
			PEEK	6	1.6	A-002014
	Small gold electrode	PEEK	3	1.6	A-002314	
Silver	Standard silver electrode		PEEK	10	5	A-002416
				6	3	A-002419
			PEEK	6	1.6	A-002011
	Small silver electrode		PEEK	3	1.6	A-002315
Palladium	Standard palladium electrode		PEEK	6	1.6	A-002019
	Small palladium electrode		PEEK	3	1.6	A-002319
Nickel	Standard nickel electrode		PEEK	6	1.5	A-002016
Copper	Standard copper electrode		PEEK	6	1.6	A-002017
			PEEK	6	3	A-012584
Iron (99.65% purity)	Standard iron electrode		PEEK	6	1.5	A-002018
			PEEK	6	3	A-012585
Carbon paste	Standard carbon paste electrode hole	depth 4 mm	PEEK	6	3	A-002210
	Small carbon paste electrode hole dep	oth 4 mm	PEEK	3	1.6	A-002223
	Carbon paste oil base 1 g				A-001010	
Boron-doped diamond	Doping level between 500 and 1000 ppm. The electrode is a 500 μ m thick disk attached to a conductive rod in brass. It is polished with an Ra < 10 nm.		PEEK	7	3	M-BDD-3

Electrode Polishing Kit

	Purpose	Catalog n°
PK-3 electrode polishing kit*		A-011975
Content		
0.05 µm polishing alumina (20 mL)	For final polishing	A-001050
1 μm polishing diamond (10 mL)	For intermediate polishing	A-002054
Glass plate (1 piece)		A-002249
Alumina polishing pad (10 pieces)	For final polishing	-
Diamond polishing pad (10 pieces)	For intermediate polishing	-
Spare parts		
Alumina polishing pad (20 pieces)	For final polishing	A-001040
Diamond polishing pad (20 pieces)	For intermediate polishing	A-001041
Emery paper UF800 (20 pieces)	For PG and PFCE electrodes	A-012611
Coarse polishing pad (20 pieces)	Rough hewn	A-001042
6 μm polishing diamond (10 mL)	For intermediate polishing	A-002053




*To refresh the working electrode surface, we recommend polishing before each measurement.

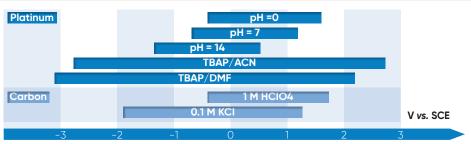

Window of electroactivity for working electrodes in different media


Counter Electrodes

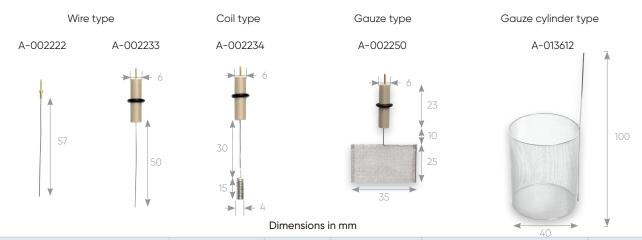

	Size/mm	Wire Ø/mm	Surf. area/cm²	Purpose	Catalog n°
Platinum (99.95 % purity)	57	0.5	~ 0.7	SVC-2, VC-4, plate material evaluating cell	A-002222
	50	0.5	~ 0.7	SVC-3	A-002233
	45 (230 coiled)	0.5	~ 3.6	RRDE, bulk electrolysis cell, SVC-3	A-002234
Gold	45 (230 coiled)	0.5	~ 3.6	RRDE, bulk electrolysis cell, SVC-3	A-012638
Nickel	45 (230 coiled)	0.5	~ 3.6	RRDE, bulk electrolysis cell, SVC-3	A-012639
Platinum gauze electrode, PEEK body 80 mesh (99.95 % purity)	25x35	0.08	~ 22.9	Bulk electrolysis cell	A-002250
Platinum gauze electrode, 54 mm wire 80 mesh (99.95 % purity)	25x35	0.08	~ 22.9	Flat cell	A-702439
Platinum electrode, 80 mesh (99.95 % purity)	40x50	0.12	~ 47.4	Bulk electrolysis cell	A-013612
Gold gauze electrode, PEEK body 100 mesh	25x35	0.07	~ 29	Bulk electrolysis cell	A-002251

Small-Size Reference Electrodes for aqueous solutions

There are two main types of reference electrodes: those suited for aqueous solutions and those for organic solutions.

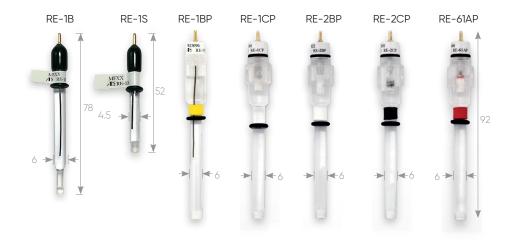

	Junction	Electrolyte	Purpose	Catalog n°
RE-1B Ag/AgCl	IPPG*	3 M NaCl	SVC-2, SVC-3, VC-4, bulk electrolysis cell, RDE/RRDE, flat cell	A-012167
RE-1S Ag/AgCl	IPPG*	3 M NaCl	SECM	A-012168
RE-1BP (Ag/AgCI)	Ceramic	3 M NaCl	SVC-2, SVC-3, VC-4, bulk electrolysis cell, RDE/RRDE, EQCM, flat cell	A-013613
RE-1CP Ag/AgCl	Ceramic	Saturated KCl	SVC-2, SVC-3, VC-4, bulk electrolysis cell, RDE/RRDE, flat cell	A-013429
RE-2BP Hg/Hg ₂ Cl ₂	Ceramic	Saturated KCl	SVC-2, SVC-3, VC-4, bulk electrolysis cell, RDE/RRDE, flat cell	A-013430
RE-2CP Hg/Hg ₂ SO ₄ , free from chloride	Ceramic	Saturated K ₂ SO ₄	SVC-2, SVC-3, VC-4, bulk electrolysis cell, RDE/RRDE, flat cell	A-013431
RE-61AP Hg/HgO main body in polyacetal resin	Ceramic	1 M NaOH	For alkaline media	A-013395
Spare parts				
Porous glass frits (CoralPor)				See p.30
Options				
RE-PV preservative vial for reference electrode, 10 mL				A-012108
Bridge tube Ø 9.0 mm (2 pieces) - Picture available next page				A-012177
Bridge tube Ø 9.0 mm (22 pieces) - Picture available next page				A-012307

^{*} Ion Permable Porous Glass


Support: How to Check Your Reference Electrode

Click or scan

BioLogic's Learning Center is a great source of information for tech-tips, theory and product information. Visit our Learning Center to find out why reference electrode maintenance is so important.

Small-Size Reference Electrodes

for non aqueous solutions

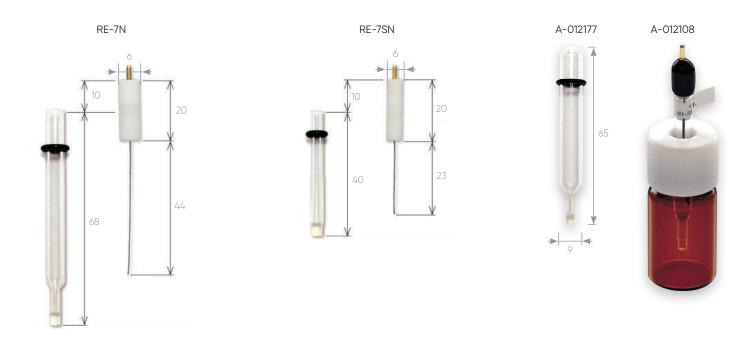

	Junction	Electrolyte	Purpose	Catalog n ^o	
RE-7N (Ag/Ag ⁺)	IPPG*	Ag/Ag ⁺ /ACN** TBAP*** (not provided)	SVC-2, SVC-3, VC-4, bulk electrolysis cell, RDE/RRDE, eQCM	A-013848	
RE-7SN (Ag/Ag ⁺)	IPPG*	Ag/Ag ⁺ /ACN** TBAP*** (not provided)	SECM	A-013849	
Spare parts					
PTFE cap with Ag wire (for RE-7N)					
Sample holder 6 mm diameter (for RE-7N) (2 pieces)					
Porous glass frits (CoralPor)				See p.30	
Options					
RE-PV preservative vial for reference electrode, 10 mL				A-012108	
Bridge tube Ø 9.0 mm (2 pieces)				A-012177	
Bridge tube Ø 9.0 mm (22 pieces)				A-012307	

- Ion Permable Porous Glass
 Acetonitrile
 Tetrabutylammonium perchlorate

Don't forget!

There is a huge amount of supporting information on <u>www.biologic.net</u>.

The BioLogic Learning Center has over 150 articles and if you want a deeper understanding, there are over 80 application notes and 50 technical notes on the field of electrochemistry alone.

biologic.net

Maintenance of Reference Electrodes

Potentials of common reference electrodes

Store your reference electrode immersed in the electrolyte

When not in use, we recommend that you keep reference electrodes in sealed, air-tight vials in order to prolong their life. The storage solution should be identical to the filling solution of the reference electrode. Preferably a cold and dark place.

Prevent contamination

To prevent contamination of the reference electrode, a bridge tube can be used.

E at 25°C
Hg/HgO / NaOH (0.1 M)
$Hg/Hg_2SO_4/K_2SO_4$ (sat)
Fc/Fc+ / TBAP (0.1 M) ACN
Ag/Ag+ / TBAP (0.1 M) ACN
Hg/Hg ₂ Cl ₂ / KCl (sat)*
Hg/Hg ₂ Cl ₂ / NaCl (sat)*
Ag/AgCl / KCl (3.5 M)
Ag/AgCl / KCl (sat)
Ag/AgCl / NaCl (sat)
NHE Normal Hydrogen Electrode * Hg/Hg ₂ Cl ₂ : Calomel

Replace the junction when needed

If you are using IPPG junctions, yellowish discoloration indicates contamination. This is caused by the absorption of organic compounds into the pores. The average pore diameter of IPPG is about 40 - 200 Å. If you are using CoralPor™ junctions, you might want to use one of the available replacement kits. The average pore diameter of CoralPor™ is about 4 - 10 nm.

	Compatible reference electrodes	Compatible bridge tubes	Content
092-VYC3	A-012167	A-012176 A-012306 A-012177 A-012307	10 glass frits (Ø 2.8 mm Coralpor™) 200 mm long heat shrink tube (Ø 3.2 mm)
092-VYC4	-	EL-C-005 EL-C-017 EL-F-004B EL-F-004 EL-A-017 EL-A-008	10 glass frits (Ø 4 mm CoralPor™) 200 mm long heat shink tube (Ø 4.8 mm)
092-VYC5	A-012168	-	10 glass frits (Ø 2.8 mm CoralPor™) 200 mm long heat shrink tube (Ø 4.8 mm)


Hydrogen Reference Electrodes

These hydrogen reference electrodes, manufactured by Gaskatel, are beneficial in that they are Hg-free and not made of glass, which extends their range of operating conditions. These reference electrodes are easy to use and robust. The hydrogen source is contained within a cartridge that is easily replaceable.

	Catalog n°
HydroFlex® Hydrogen Reference Electrode	G-HYDROFLEX
HydroFlex® Starter Kit*	G-HYDROFLEX-KIT
MiniHydroFlex® Hydrogen Reference Electrode	G-MINIHYDROFLEX
Options	
HydroFlex® Hydrogen Cartridge (x4)	G-HYDROFLEX-CARTRIDGE
MiniHydroFlex® Hydrogen Cell (x2)	G-MINIHYDRO-CELL
Bridge tube for use of HydroFlex®with EL-ELECTRO	EL-A-017
Bridge tube for use of HydroFlex®with EL-CORR	EL-C-017
Bridge tube for use of HydroFlex®with EL-FLAT	EL-F-004B

	Operating conditions
Cartridge lifetime/months	12 (MiniHydroFlex®), 6 (HydroFlex®)
Temperature/°C	-20 to 120 (PTFE body)
pH range	-2 to 16

The HydroFlex® electrode is compatible with BioLogic's electrochemical cells. The MiniHydroFlex® is dedicated to the FlexCell®.

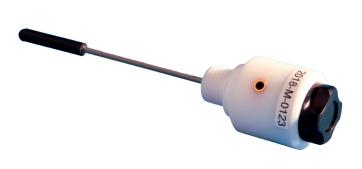

Part	Materials
Body	PTFE, PP
Shaft	PEEK, PTFE
Cap and cartridge	PC or PVC
Measuring electrode	Pt, Pd
PP	Polypropylene
PTFE	Polytetrafluoroethylene
PVC	Polyvinyl chloride
PC	Polycarbonate
PEEK	Polyether ether ketone
Pt	Platinum
Pd	Palladium



Dimensions in mm

HydroFlex® Starter Kit Content
1x HydroFlex® hydrogen reference electrode (incl. 1 hydrogen cartridge)
1x connector lead
1x hydrogen cartridge
1x cartridge wrench
1x operating time wrench



Connection Accessories.

High-Temperature Extension Cables

These extension cables are intended to be connected between the instrument's cell cable and the cell.

Available in two different sizes, 1.3 m and 2.5 m long, and with a maximum temperature of 150°C, these cables are compatible with:

- SP-50e, SP-150e, VSP, VSP-3e, VMP-3e and boosters for Essential range instruments
- SP-200, SP-240, SP-300, VSP-300, VMP-300, and boosters for Premium range instruments
- BCS-905, BCS-910 and BCS-915

	1.3 m	2.5 m
Temperature /°C	-40 to 150	
Max current /A (2 mm connectors)	2	
Max current /A (4 mm connectors)	20	
Cable diameter /mm	12.7	
Catalog n°	092-25/1	092-25/2



Glove Box Cables



Hermetic cell cable for glove boxes

As standard, the potentiostat and the booster are provided with a 1.5 m long cell cable. The cable connected from the booster to the potentiostat is 0.8 m long for VMP-3e based instruments.

For some applications, the user may need different cable lengths. For this reason, longer cables are available (for more information, contact your sales representative).

For applications carried out in glove boxes, special glove box cell cables are also available.

Setup Connection

Bad connections can affect measurements (stability of potentiostat, artefacts etc).


In order to optimise your setup, we recommend you use the accessories described in this section.

	Essential	Premium
Catalog n°	092-23/5	094-101/6 (standard cable)
		094-101/8 (low current cable)
Content		
Feedthrough type	12 pins	25 pins*
Inside glovebox (length: 1 m)	Cable with 2 mm connectors on one side and 12-pin Jaeger connector on the other side	Cable with electrometer on one side and 25-pin connector on the other side
Outside glovebox (length: 1.5 m)	Cable with Sub-D 25 connector on one side and 12 pin Jaeger connector on the other side (length 1.5 m)	Cable with Sub-D 25 connector on one side and 12 pin Jaeger connector on the other side (length 1.5 m)
Requirement		
Hole diameter needed in glovebox /mm	27	45

^{*}Two feedthrough seals : one installed in the glove box wall, the other dedicated to the channel board.

Multi-Electrode Investigation Cables

For the Essential product range, we offer several options to facilitate the use of connection cables when multi-electrode experiments are performed: for example, RRDE experiments or corrosion experiments on several samples using the same reference and the same counter electrode.

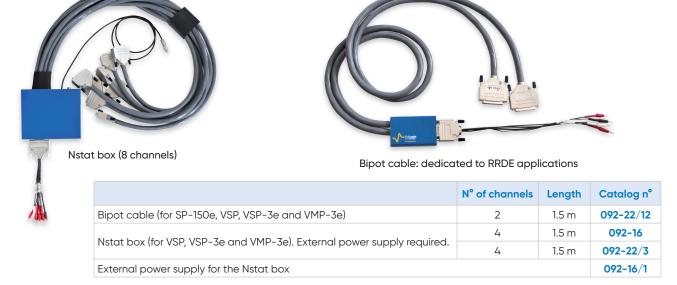

Connectors

	Picture	Size	Colors	Catalog n°
Aligator clips	7 7	2 mm	8x Red, 8x Blue, 8x White, 4x Black	092-1001/40
	7. 3.	4 mm	3x Red, 3x Blue, 3x Black	092-1001/41
Receptacles		2 mm	20x Red, 20x Blue, 20x White, 15x Black	092-1001/42
		4 mm	12x Red, 12x Blue, 12x Black	092-1001/43
Banana plugs		2 mm	10x Red, 10x Blue, 10x White, 10x Black	092-1001/44
		4 mm	8x Red, 8x Blue, 8x White	092-1001/45
Adaptors	1111	2 mm receptacle to 4 mm plug adapter	6x Red, 6x Blue, 4x White, 4x Black	092-1001/46
		2 mm receptacles to 4 mm banana plugs	16x Red, 16x Blue	092-1001/47
		4 mm receptacles to 2 mm banana plugs	20x Black	092-1001/48

Connection Kits

	Content	Colors	Catalog n°
For standard board	4 alligator clips of 2 mm 3 receptacles of 2 mm	red, blue, white, black red, blue, white	092-1001/30
For booster board	3 alligator clips of 2 mm 2 alligator clips of 4 mm 3 receptacles of 2 mm 2 receptacles of 4 mm	red, blue, white red, black red, blue, white blue, white	092-1001/31
For HCV-3048, FlexP 0060, 0160, CC4-60A and CC8	Contains: 2x6 mm receptacles 4 lugs with 4 mm receptacles	red, blue red, blue	094-110/CNT
For FlexP0012 and CC4-200A	Contains: 2 Amphenol 8 mm receptacles, 4 lugs with 4 mm receptacles	red, black red, blue	093-200/CNT





093-200/CNT

External Device Connection





Catalog n° DB9-8 BNC connector for auxiliary I/O 092-22/1 IS1 isolation module for auxiliary I/O for Premium based 094-081/5 PT100 temperature probe, to be connected to the 092-22/13 auxiliary I/O, temperature range: -50°C to 250°C, Dimensions: 3 x 20 mm, Length of cable: 2.5 m, Accuracy: ±1°C For temperature measurement in air PT100 temperature probe for temperature EL-C-014 measurement in solution with Sub-D 9 connector PT100 temperature probe for temperature EL-C-014/1 measurement in solution with triad connector

Test Boxes

	Description	Catalog n°
Test Box 2	Several circuits with high precision resistors for calibration and validation	092-22/6
Test Box 3	Three circuits: linear, two non-linear systems (Tafel & passivating) for teaching and demonstration	092-22/7



Test Box 2 Test Box 3

Faraday Cage


To avoid any external perturbations, especially for low current applications, we recommend using a Faraday cage.

Please note that in order to ensure that the cage is fully functional, it must be earthed by connecting it to the ground (this is done via a green plug on the instrument's rear panel).

	Catalog n°
FC-45 Faraday cage, 450x450x450 mm	094-084/1
Cell stand for FC-45	094-084/2



Spectroelectrochemistry.

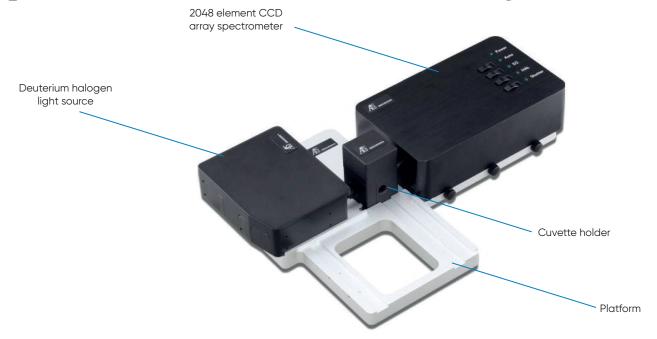


Spectrometer System

Spectroelectrochemistry (SEC) can be useful to elucidate electrochemical reaction mechanisms. The spectroelectrochemical kit is made up of three parts (spectrometer, light source and cuvette holder).

The spectrometer is equipped with a trigger to synchronise electrochemical and spectroscopic measurements.

medadiemento.	
	Information
Detector	2048 element linear silicon CCD array
Full description	SEC2021-025-DUVN
Detector range/nm	200 - 1025
Grating	Blaze wavelength (300 nm)
Slit/µm	25
Wavelength resolution/nm	1.3
Fiber connector	SMA905, core diameter: 600 µm NA=0.22
Interface	USB2.0
Operating system	Windows 10/11 (32bit/64bit)
Dimensions (HxWxD)/mm	32x86x110

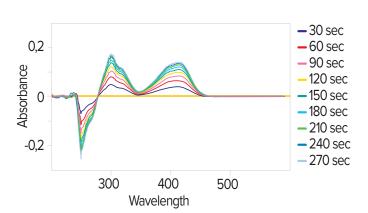

	Catal	og n°
SEC2020 spectrometer kit	A-013609)
Content		
SEC2021 Spectrometer (x 1)		
SEC2022 Deuterium halogen light source (x 1		
SEC2023 Cuvette holder (x 1)		
SEC2024 Platform (x 1)		
AC adaptor (x 1)		
Power cable (x 1)		
USB cable (x 1)		
Collimator (x 2)		
Fibre collimator (x 1)		
Platform screw (x 7)		
External device connection trigger cable (x 1		
Light source control trigger cable (x 1)		
Plastic cuvette (x 1)		
SMA905 adaptor for light source (x 3)		
SMA905 adaptor for light shielding (x 2)		
Hexagon wrench 0.89 mm (x 1)		
Hexagon wrench 1.50 mm (x 1)		
Software (USB memory) (x 1)		
Waterproof box (x 1)		
Quick manual, wavelength calibration data sheet and warranty certificate are also inclu-		est data
Option		
Connecting cable to synchronise the SEC202 with BioLogic instrument.	092-22/1	1

The SEC2020 spectrometer system uses the Czerny-Turner optical mount. This system is a M-shaped structure symmetric to the grating (4) and is an optical system with extremely small aberration.

Light source structure




	Information
Light type	Deuterium halogen light source
Wavelength range/nm	200 - 1700
Stability	<0.1%
Drift/h	0.25%
Bulb life/h	>1000 (D2 lamp) >2000 (halogen lamp)
Fiber connector	SMA905
Size (HxW×D)/mm	46x100×165


- 1. SMA905 Connector
- 2. Slit
- 3. Collimating mirror
- 4. Grating
- 5. Focus mirror
- 6. 2048 element CCD array

Utilization modes


	Information
Transmittance	'
Absorbance/Transmittance	- Concentration of chemicals (solution) - Polymer extrusion processes - DNA quantification
Reflectance	- Freshness testing - Film thickness/composition (quality control) - Activation energy of photocatalytic species - Textile quality control
Fluorescence	
Fluorescence	- Marine organisms - Biology (DNA, protein, cell proliferation assay, histamine-analysis, algae monitoring) - Environmental fields (waste water analysis, ground water trace studies, hydrocarbon detection, dissolved oxygen) - Plant efficiency (plant physiology, plant breeding, horticulture, agronomy, agrochemicals, pollution studies, forestry, ecology) - Tissue diagnosis
Scattering	Oil concentrations of oil/water systems Raman spectroscopy Physical transition phenomena (e.g: melting point, glass transition crystallization temperature)
Irradiance	
Emission	Astronomy (e.g., spectra of Hale-Bopp, plasma monitoring) In situ metal monitoring Luminescence (Photoluminescence, Electroluminescence), LED & laser wavelength

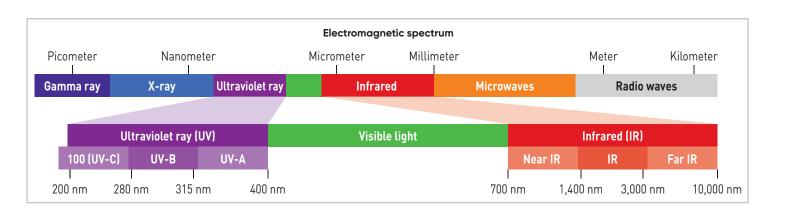



Spectroelectrochemical Cell

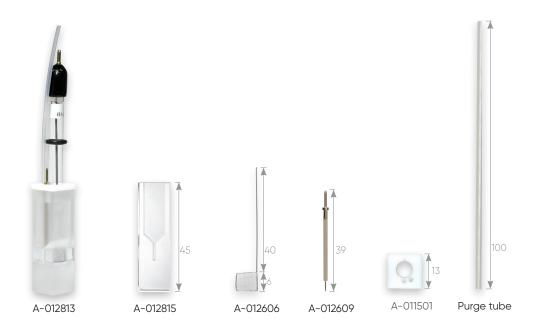

		Platinum	Gold	
0.5 mm	Kit	A-012813	A-012814	
	Content			
	Platinum counter electrode	A-012609		
	Thin layer quartz glass cell	A-012815		
	PTFE cap	A-	011501	
	Purge tube (ETFE, 100 mm)	A-012606 (80 mesh, height 5 mm) - (100 mesh, height 5 mm)		
	Gauze working electrode			
1 mm	Kit	A-013510	A-013511	
	Content			
	Platinum counter electrode	Α-	012906	
	Thin layer quartz glass cell	Α-	012907	
	PTFE cap	A-	011501	
	Purge tube (ETFE, 100 mm)		-	
	Gauze working electrode	A-011498 (80 mesh, height 6 mm)	A-012017 (100 mesh, height 6 mm)	
Options				
RE-1BP A	g/AgCl reference electrode	A-	013613	
RE-7N no	n-aqueous reference electrode	A-013848		
Purge tuk	pe (ETFE), 1 m	A-010537		

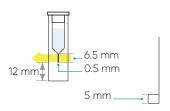

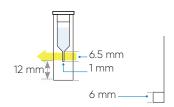

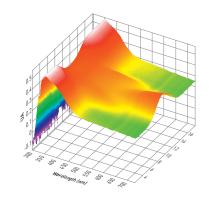









For more information visit biologic.net



Application notes



Learning center



Tutorials



Videos



We serve our customers worldwide through our subsidiary offices and our extensive distribution network.

Headquarters BioLogic SAS

4, rue Vaucanson 38 170 Seyssinet-Pariset France

Phone: +33 476 98 68 31 Fax: +33 476 98 69 09



Subsidiaries offices BioLogic USA, LLC

USA

Phone: +1 865 769 3800

BioLogic Science Instruments GmbH

Germany

Phone: +49 (0)551 38266900

BioLogic Science Instruments Ltd

United Kingdom

Phone: +44 (0) 333 012 4056

BioLogic Spain

Spain

Phone: +34 681 357 873

BioLogic Pvt Ltd

India

Phone: +91 022 46055588

BioLogic Singapore

Singapore

Phone: +65 92335838

Always by your side. Wherever you are.

We here at BioLogic pride ourselves on the quality and robustness of our instruments. However if you, for whatever reason, encounter a problem with your instrument or accessories, our global support network will help find you a solution quickly and effectively.

If you need more information, or perhaps a little inspiration, you can browse our ever-growing support database with hundreds of Learning Center articles, application/technical notes and support videos at www.biologic.net.

